


## Standard Operating Procedure (SOP) for PV360.2.0.pl - Example Image Stitching (mouse full body)

See SOP Study Initiation Protocol for more information on how to start image acquisition.


1. **Imaging:** Drag to instruction window the scan of interest. In this example, we'll work with the Multislice Localizer

- **Geometry Card:** we can adjust slices (number, orientation, thickness, gap, distance) and bed positions (should be >1 for stitching)
- More details regarding the position and bed shift distances can be visualized if clicking on the "..."

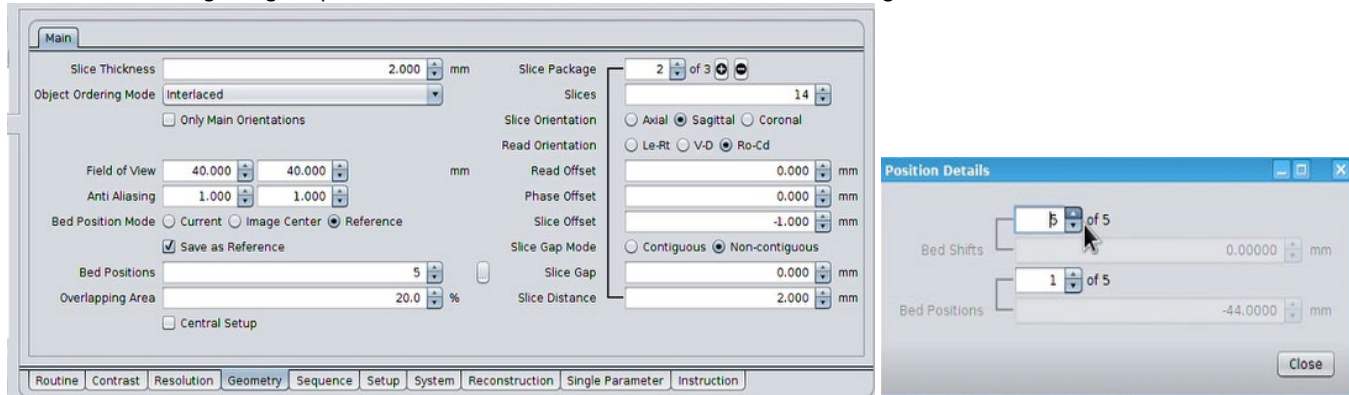


2. **Stitching:**


- After acquisition: right click in Palette>Explorer>Scan EX/Data Reconstruction (P1) → Create Image Series... > MRI Data Reconstruction > OK

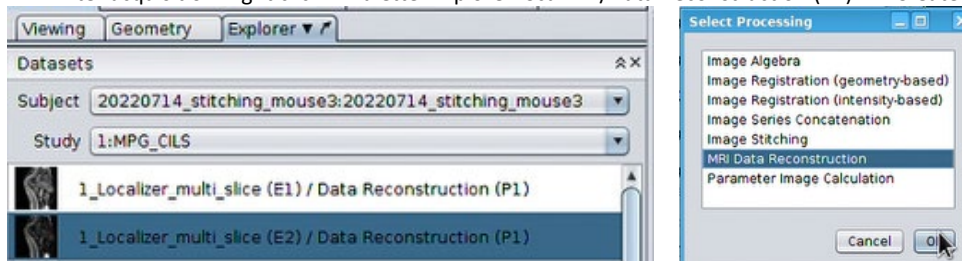


- Repeat in the resulting scan ("MRI Data Reconstruction (P2)") but this time select "Image Stitching"

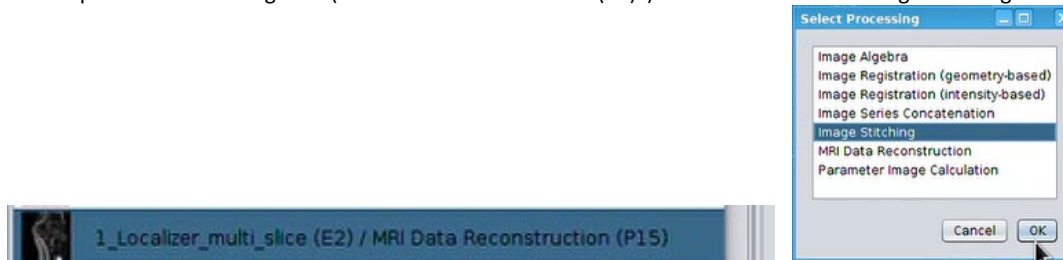

- The stitching reconstruction menu will appear (these are the suggested parameters), select Execute when done

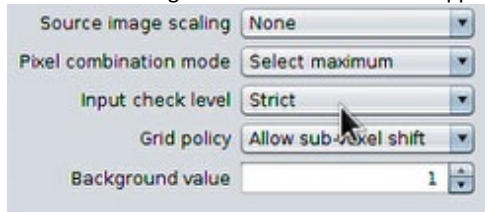

- Check the images to make sure the resulting stitching makes sense before saving them

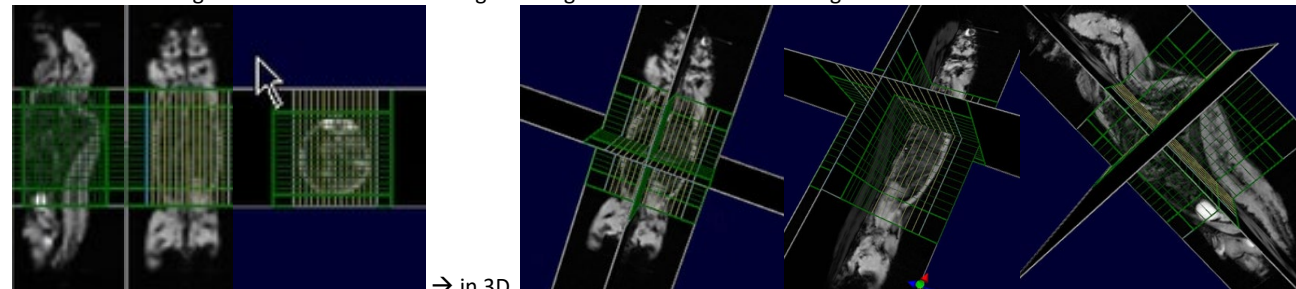

3. **Back-up & Export of image data** (automatically in the large storage drive)

- Additional options to save as raw Bruker files or DICOM or Nifti via **Palette > Explorer > right click>'Export to...'** 'Topspin'/'Dicom'/'Nifti'....